

Specification

Long Pulsed Nd:YAG Allux Nd:YAG

Wavelength	1064nm
Laser Type	Long Pulsed Nd:YAG
Beam Delivery	Lens-coupled optical fiber with Handpiece
Aiming Beam	5 mW @532 nm
Energy	90J
Fluence Max.	1273J/cm ²
Pulse Width (ms)	0.5ms-100ms
Pulse Rate (Hz)	0.5-7Hz
Display System	10.2 " Touch LCD GUI
Spot Size	3, 5, 7, 10, 12, 15, 18, 20, 22 mm
Emission Control	Footswitch
Cooling Type	Air Cooling / Gas Cooling
Dimensions (mm)	1000mm (H), 910 mm (D), 380 mm (W)
Weight (kg)	88 kg
Electrical rating	10 kVA (AC 220 ~ 230 V,50/60 Hz)

Alexandrite Allux Alex

Wavelength	755nm
Laser Type	Alexandrite
Beam Delivery	Lens-coupled optical fiber with Handpiece
Aiming Beam	5 mW @532 nm
Energy	60J
Fluence Max.	849J/cm ²
Pulse Width (ms)	0.5 ms-100ms
Pulse Rate (Hz)	0.5-7Hz
Display System	10.2 " TouchLCD GUI
Spot Size	3, 5, 7, 10, 12, 15, 18, 20, 22mm
Emission Control	Footswitch
Cooling Type	Air Cooling / Gas Cooling
Dimensions (mm)	1000mm (H), 910 mm (D), 380 mm (W)
Weight (kg)	88 kg
Electrical rating	10 kVA (AC 220 ~ 230 V, 50/60 Hz)



RM, 1412, SJ Techno ville, 278, Beotkkot-Ro, Geumcheon-Gu, Seoul, Korea
Tel: 82-2-3281-9091 Fax: 82-2-3281-8011 E-mail: amimed@amimed.co.kr



Allux Dual


Alexandrite & Long Pulsed Nd:YAG Laser

Alexandrite & Long Pulsed Nd:YAG Laser

World best laser for Hair removal



Before & After

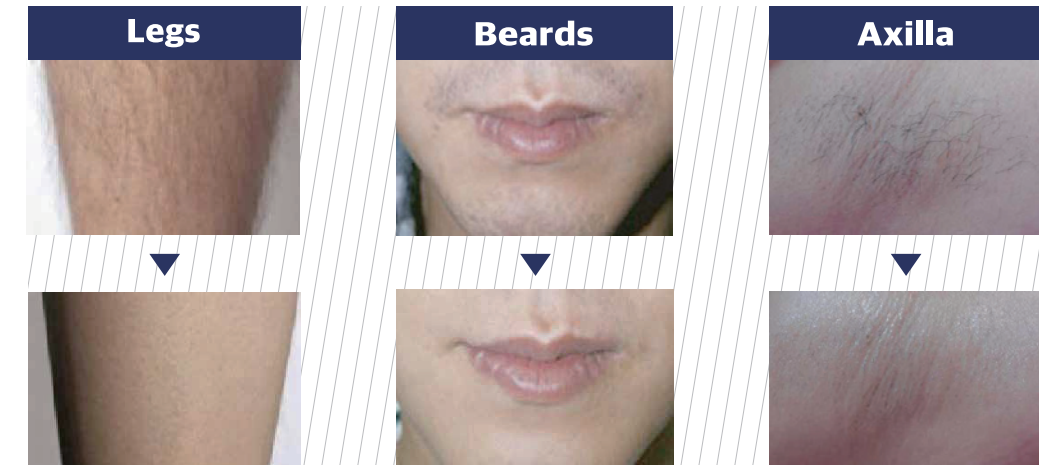

Strong points

Laser beam is transferred through fiber to the wanted area precisely without energy loss. So, it is very convenient to use and the treatment time can be quite shortend. Cryo Air Cooling system is installed. So, it can treat patients without pain and also, downtime can be sharply shortend.

Features

- 01 The gas capacity for one cylinder is very big (900g) comparing other products. With one full filled cylinder, 25,000 shots possible in 4ms condition. Also, it can be refilled.
- 02 Optimized automatic temperature controlling system. Temperature of cooling water can be adjusted automatically.
- 03 Subdivided spot sizes can make the treatment very convenient and precise.
- 04 Easy replacement of water and filter with one-touch system available.



• Hand piece

Removable handpiece



Zoom handpiece

